

VOCO Profluorid® Varnish

YOUR CHOICE

Indications

- Treatment of hypersensitive teeth
- Sealing of the dentinal tubules on sensitive root surfaces or following cavity preparation

Advantages

- Easy to apply on wet surfaces
- Pleasant melon flavour
- Aesthetic, white-transparent tooth colour
- Convenient SingleDose, product and applicator in hygienic disposable packaging
- Fast desensitisation and fluoride release (5% NaF \approx 22,600 ppm fluoride)



Presentation

REF 1268 50 SingleDose 0.25 ml, accessories
REF 1269 50 SingleDose 0.40 ml, accessories



Presentation

REF 1267 Tube 10 ml



VOCO GmbH
P.O. Box 767
27457 Cuxhaven
Germany
Tel.: +49 (0) 4721-719-0
Fax: +49 (0) 4721-719-140
info@voco.de
www.voco.com

Available from:

VOCO Profluorid® Varnish

DESENSITISING VARNISH WITH FLUORIDES

VOCO Profluorid® Varnish

HIGHLY EFFECTIVE

Hypersensitivity

Sensitive teeth are a widespread problem for a large number of people. Usually, the condition causes only transient pain which patients do not always report to their dentists.


Hypersensitivity is mainly caused by exposed dentinal tubules. Therefore, the aim of treatment is to seal these tubules.

Indications

VOCO Profluorid® Varnish is a fluoride-containing varnish for treating hypersensitivity. Dentinal tubules exposed on the necks of teeth or after cavity preparation are effectively sealed.


Components and effect of VOCO Profluorid® Varnish

VOCO Profluorid® Varnish is a colophony-based varnish containing 5% sodium fluoride (22,600 ppm fluoride). The fluoride ion, together with the calcium ions accumulated in the tubules, causes a precipitation of calcium fluoride, which effectively seals the tubules. According to a user survey, over 85% of patients were completely free of pain after only one treatment with **VOCO Profluorid® Varnish**.

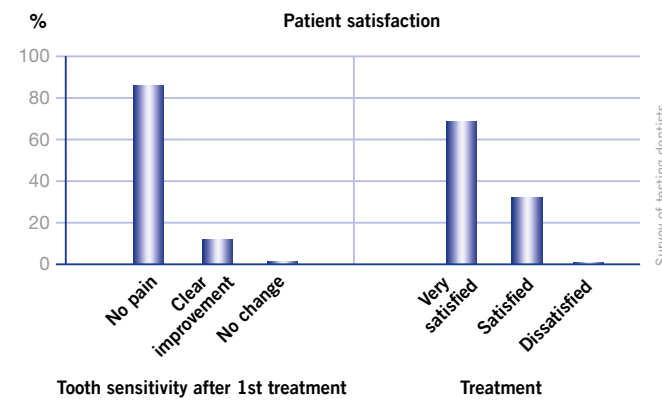

As well as rapidly sealing the dentinal tubules, **VOCO Profluorid® Varnish** also causes calcium fluoride deposits to form on the tooth surface. These deposits protect the tooth from acid attacks, promote remineralisation and contribute to the long-term formation of fluoroapatite.

Xylitol, which has a proven cariostatic effect in addition to its taste-enhancing properties, is also added to the varnish.

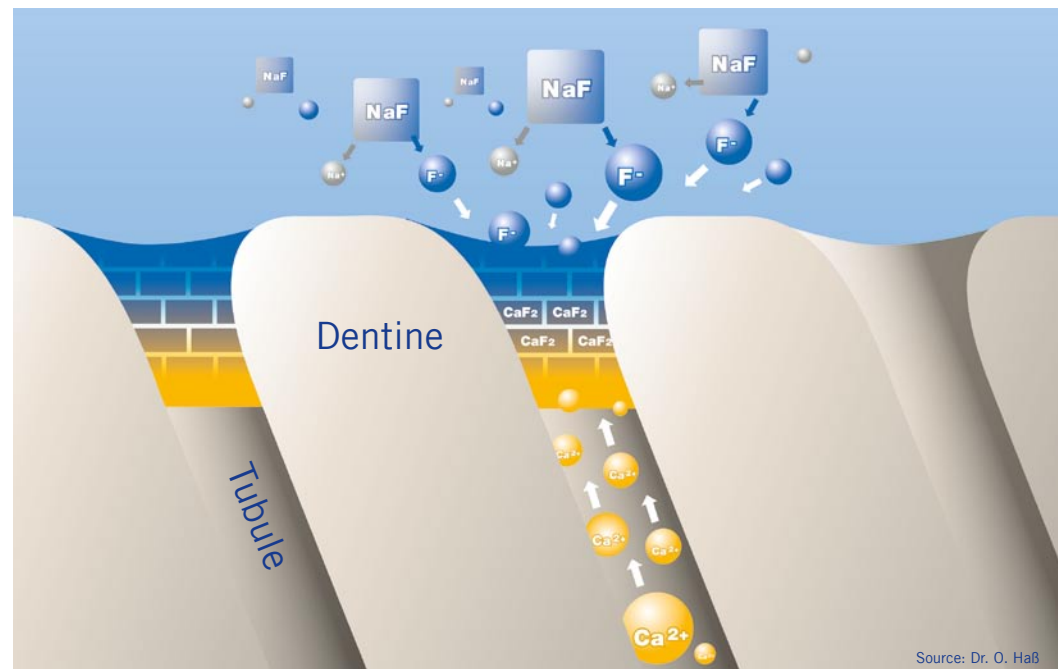

Source: Dr. O. Haß

VOCO Profluorid® Varnish





SIMPLE AND AESTHETIC

Easy to apply

Thanks to its colophony matrix, **VOCO Profluorid® Varnish** adheres well even to wet surfaces. It is not necessary to dry the tooth surfaces. This improves cooperation of patients with low levels of compliance (e.g. children).

Aesthetics and taste

VOCO Profluorid® Varnish is a tooth-shaded varnish which does not impair aesthetics. Melon flavour is added to the varnish in addition to xylitol to improve taste and aroma. This considerably increases acceptance, especially with children.

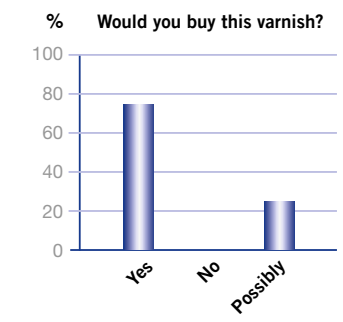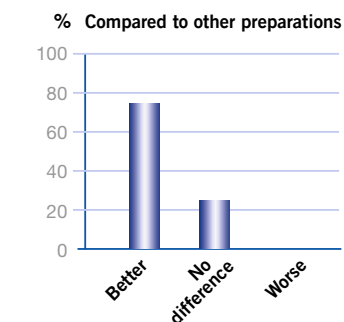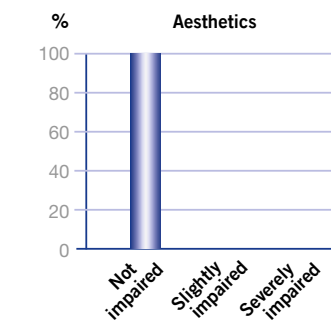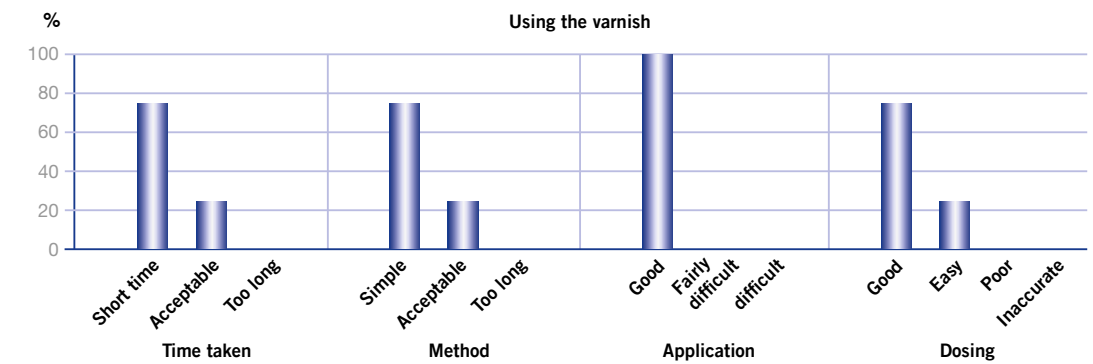

Opinion of dentists

Dentists using the product all confirmed its good handling and aesthetic properties. A vast majority also reported positive differences between Profluorid Varnish and competing products. **VOCO Profluorid® Varnish** was equally well accepted by the patients. After the treatment, nearly all of the patients reported that they were either symptom-free or had a significantly reduced level of sensitivity.

VOCO Profluorid® Varnish: much more than just good taste!

The SingleDose advantages

VOCO Profluorid® Varnish SingleDose comes in a blister pack containing a brush and the varnish. This enables the product to be applied under perfectly hygienic conditions. The quantity of **VOCO Profluorid® Varnish** SingleDose is pre-dosed in the blister. The 0.25 ml version is adapted to the needs of a deciduous dentition, the 0.4 ml blister allows to treat an adolescent or adult dentition. This reduces the risk of over-dosing.

Simple and aesthetic desensitising with VOCO Profluorid® Varnish:

